


ENDOSCOPY EICKVIEW LED VIDEOENDOSCOPES


Versatile & innovative

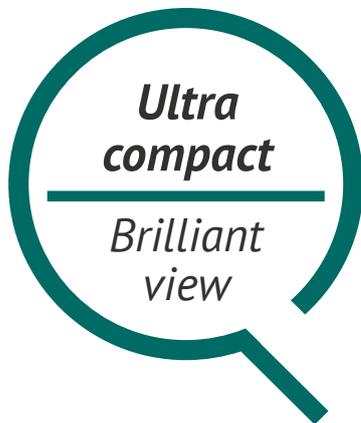

Endoscopy

VERSATILE & INNOVATIVE EICKVIEW LED VIDEOENDOSCOPES

Due to their extremely compact dimensions, the EickView LED Video Endoscopes are the perfect solution for broncho-, gastro-, rhino-, tracheo-, uretero- and vaginoscopy in dogs and cats and for the examination of the upper respiratory tract in horses.

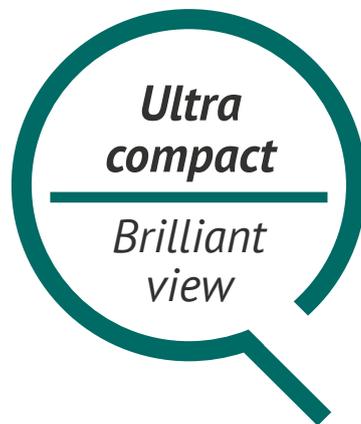

High mobility

Stand-alone system including powerful LED light source and high resolution camera.

Easy documentation

Images and videos can be recorded and archived on the integrated SD card.

Perfect image quality

New CMOS Chip technology provides high-resolution images and videos.

Many possibilities

The integrated working channel allows taking biopsies, extracting foreign bodies, flushing etc.

Simple cleaning

Simple and clear construction allows accurate cleaning in no time.

“We use our EickView LED Videoendoscope regularly, it is an extremely versatile instrument that provides valuable and high quality data for our practice.”

revolutionary



ENDOSCOPY LED VIDEOENDOSCOPES

EickView 60S

- ▶ Working Length: 60 cm
- ▶ Outer Ø: 3.3 mm
- ▶ Working channel Ø: 1.2 mm
- ▶ Deflection: 180° up, 130° down

306009

EickView 60

- ▶ Working Length: 60 cm
- ▶ Outer Ø: 4.2 mm
- ▶ Working channel Ø: 1.3 mm
- ▶ Deflection: 180° up, 130° down

306011

EickView 70

- ▶ Working Length: 70 cm
- ▶ Outer Ø: 6.0 mm
- ▶ Working channel Ø: 2.0 mm
- ▶ Deflection: 180° up, 130° down

306010

EickView 100E

- ▶ Working Length: 100 cm
- ▶ Outer Ø: 8.0 mm
- ▶ Working channel Ø: 2.8 mm
- ▶ Deflection: 150° up, 130° down

306166

EickView 150E

- ▶ Working Length: 150 cm
- ▶ Outer Ø: 8.0 mm
- ▶ Working channel Ø: 2.8 mm
- ▶ Deflection: 130° up, 130° down

306164



306009

-
- ▶ Folding 3.5" monitor for video and picture playback
 - ▶ Recording and archiving on integrated SD card
 - ▶ Can be connected to a larger monitor
 - ▶ File formats: JPG, AVI
 - ▶ Integrated powerful LED light source
 - ▶ Rechargeable Battery with up to 3 hours of operation
 - ▶ USB port and video out
 - ▶ Depth of field: 5 – 80 mm (EickView 60S: 5 – 50 mm)
 - ▶ Direction of view: 0°
 - ▶ Field of view: 110°

306009 – 306166

Endoscopy

ACCESSORIES FOR EICKVIEW LED VIDEOENDOSCOPES



BIOPSY FORCEPS FOR EICKVIEW ENDOSCOPES

- 306311 For EickView 60S/60
- 306314 For EickView 70
- 306319 For EickView 100E
- 306316 For EickView 150E



ALLIGATOR FORCEPS FOR EICKVIEW ENDOSCOPES

- 306211 For EickView 60S/60
- 306214 For EickView 70
- 306219 For EickView 100E
- 306216 For EickView 150E



RAT TOOTH FORCEPS FOR EICKVIEW ENDOSCOPES

- 306431 For EickView 60S/60
- 306435 For EickView 70
- 306430 For EickView 100E
- 306434 For EickView 150E



GRASPING FORCEPS FOR EICKVIEW ENDOSCOPES

- 306340 For EickView 60S/60
- 306346 For EickView 70
- 306341 For EickView 100E
- 306342 For EickView 150E



BASKET RETRIEVAL DEVICE FOR EICKVIEW ENDOSCOPES

- 306336 For EickView 60S/60
- 306332 For EickView 70
- 306338 For EickView 100E/150E



POLYPECTOMY SNARE FOR EICKVIEW ENDOSCOPES

- 306426 For EickView 70
- 306421 For EickView 100E
- 306423 For EickView 150E



CYTOLOGY BRUSHES FOR EICKVIEW ENDOSCOPES

- 306415 For EickView 60S/60
- 306416 For EickView 70/100E
- 306418 For EickView 150E



CLEANING BRUSHES FOR EICKVIEW ENDOSCOPES

- 306444 For EickView 60S/60
- 306440 For EickView 70
- 306445 For EickView 100E
- 306442 For EickView 150E



SILICONE TUBES FOR EICKVIEW ENDOSCOPES

- 306353 For EickView 60
- 306354 For EickView 100E
- 306352 For EickView 150E